


Sumps

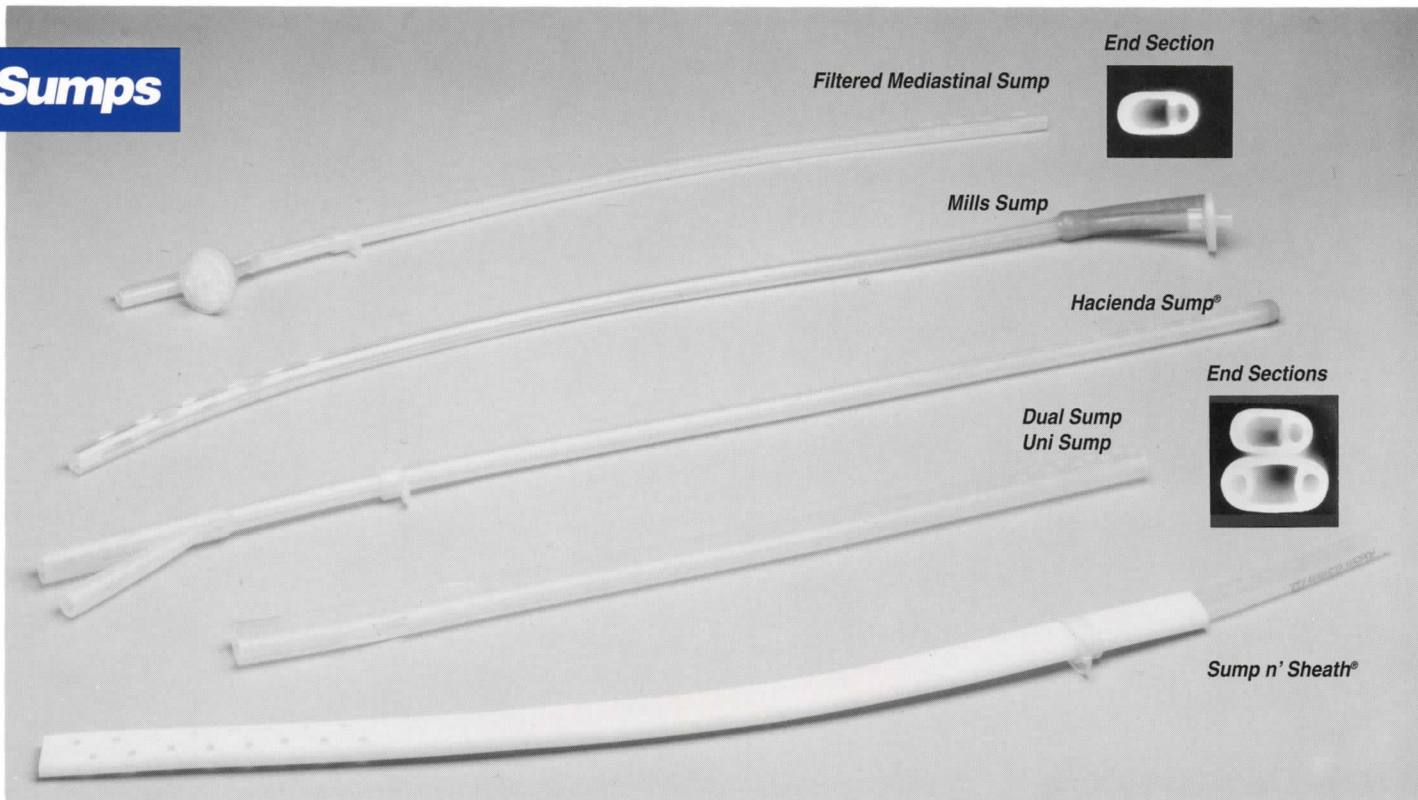

Axiom® Sumps

Axiom® silicone sump drains have the largest drainage capacity of any on the market. The multiple sump lumens create high internal flow rates that accelerate fluid removal without applying excessive or traumatic suction to delicate tissues.

The drainage eyes are of optimum size and are spaced over the maximal length to ensure evacuation of the entire drainage bed.

Filtered Mediastinal Sump

- Sump lumens create high flow rates to accelerate fluid removal
- Does not permit pressure build-up in the mediastinum

Cat. No. w/o Clot Stop®	Cat. No. With Clot Stop®	
111188	111188AT	Filtered Mediastinal Sump, 20", 8.5mm I.D., Large Eye
111189	111189AT	Filtered Mediastinal Sump, 18", 6.5mm I.D., Large Eye

Mills Filtered Mediastinal Sump

- Medication port with capped filter
- Round configuration

111167	111167AT	Mills Mediastinal Sump, 18FR
111169	111169AT	Mills Mediastinal Sump, 28FR

Hacienda Sump®

- Superior alternative to round PVC drains
- Available in pediatric and adult sizes

111162	111162AT	Hacienda Sump®, 18FR
111166	111166AT	Hacienda Sump®, 28FR

Interpleural Anesthesia Catheter

- Routine chest tube with large primary lumen, depth markings, and capped tip
- Smaller integral lumen for application of fluids (i.e. pain management)
- Three way stop-cock allows multiple options while preventing air entry

111418	111418AT	18 FR Interpleural Anesthesia Catheter - 2 lumens plus 3-way stop cock
112028	112028AT	28 FR Interpleural Anesthesia Catheter - 2 lumens plus 3-way stop cock

Abdominal Dual and Uni Sumps

- Designed with special eye patterns to eliminate the possibility of entrapped tissue
- Large internal diameter and high internal flow rates provide greatest drainage capacity of any abdominal sump available
- Available with suture rings

Cat. No. w/o Clot Stop®	Cat. No. With Clot Stop®	
111176	111176AT	Uni-Sump, 18", 6.5mm I.D., Small Eye
111177	111177AT	Uni-Sump, 12", 6.5mm I.D., Small Eye
111178	111178AT	Uni-Sump, 18", 6.5mm I.D., Large Eye
111186	111186AT	Filtered Dual-Sump, 20", 8.5mm I.D., Small Eye
111187	111187AT	Filtered Dual-Sump, 14", 8.5mm I.D., Small Eye
111196	111196AT	Dual-Sump, 20", 8.5mm I.D., Small Eye
111197	111197AT	Dual-Sump, 14", 8.5mm I.D., Small Eye
111198	111198AT	Dual-Sump, 20", 8.5mm I.D., Large Eye

Sump n' Sheath®

- Innovative drain consisting of a high capacity sump enclosed in a soft silicone sheath
- Sheath creates a pocket of negative pressure causing secretions to be swept into drain
- Air flow between sheath and drain has a cushioning effect, isolating tissue from damaging suction and reducing contact pressure

111155	111155AT	Sump n' Sheath®, 9mm - 14"
111156	111156AT	Sump n' Sheath®, 7mm - 18"
111157	111157AT	Sump n' Sheath®, 9mm - 20"

Sump Filter

111010A	Sump Filter, 99.96 D.O.P Rating
---------	---------------------------------

Suture Rings

777070	Uni-Sump Suture Ring
999090	Dual-Sump Suture Ring

All sumps and sump filters are packaged 10/box.