


SCANLAN®  
SURG-I-LOOP®  
SILICONE LOOPS

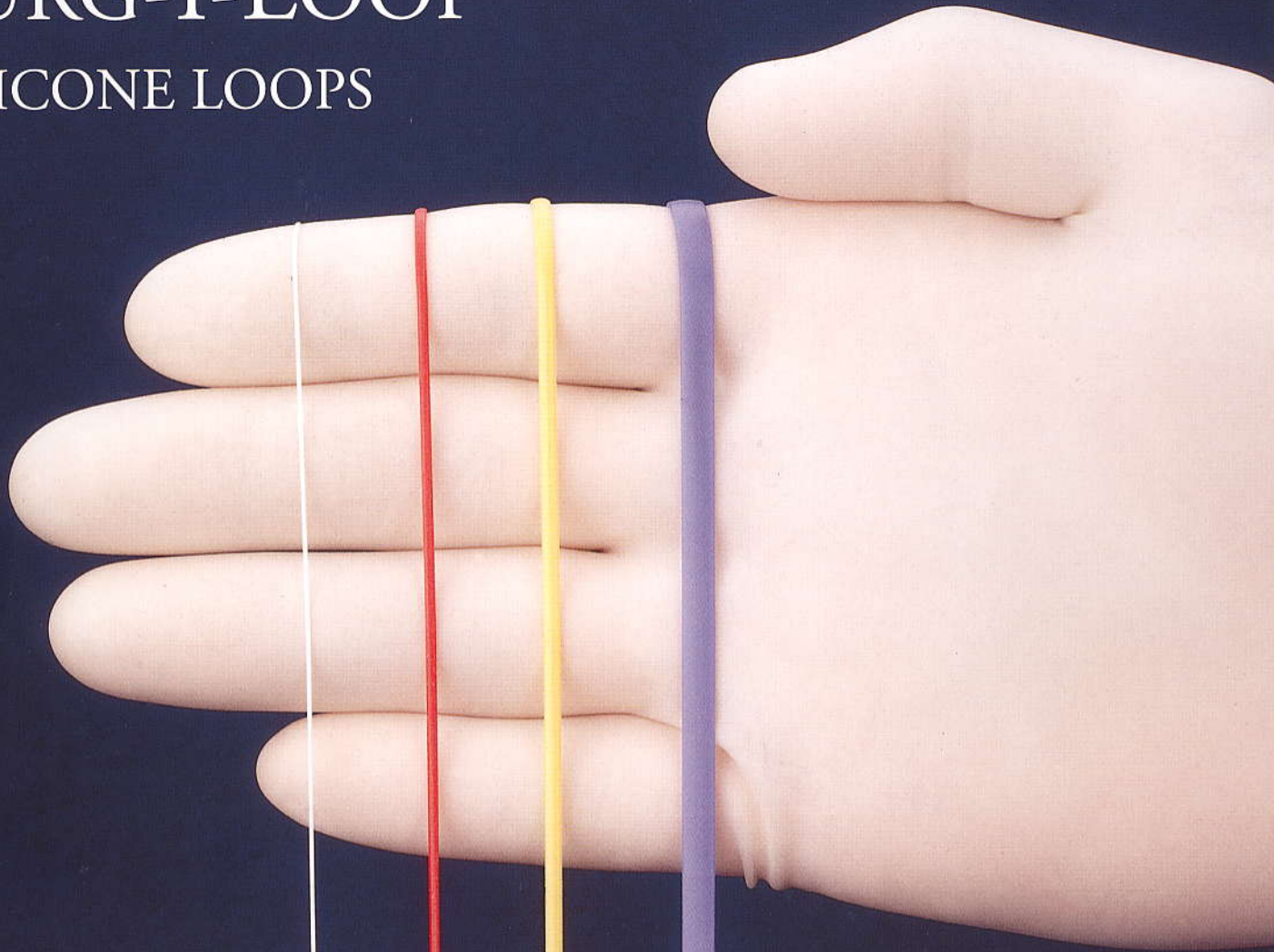

SCANLAN ★

INTERNATIONAL

*Surgical products designed and manufactured by the Scanlan family since 1921*

# SCANLAN® Surg-I-Loop® Silicone Loops

Surg-I-Loop® silicone loops are designed to aid the surgeon in complicated surgical fields by providing retraction, occlusion and identification of arteries, veins, nerves and ureters.


Manufactured of high quality, medical grade silicone to retain preset tension, replacing less appropriate retraction materials such as rubber bands, tubing, or umbilical tape.

## SCANLAN® Surg-I-Loop® Silicone Loops are...

- Manufactured from high quality, medical grade, silicone
- Designed with an oval cross-section to permit maximum surface-to-surface ratios (minimizing the trauma to vessels)
- Designed with a smooth exterior surface to slide gently over vessels
- Manufactured to retain preset tension
- Non-porous
- Sterile packaged
- Single-use
- Radiopaque
- 100% Latex Free

Surg-I-Loop® silicone loops are available in four different colors. The various colors allow easy recognition and less time consuming identification of elements during surgery. For example, red may be used for arteries, blue for veins, yellow for ureters, and white for nerves and tendons.

Surg-I-Loop® silicone loops are available in micro, mini, maxi, and super maxi sizes for use in all types of surgeries including the thoracic, neuro, plastic, orthopedic, vascular specialties.

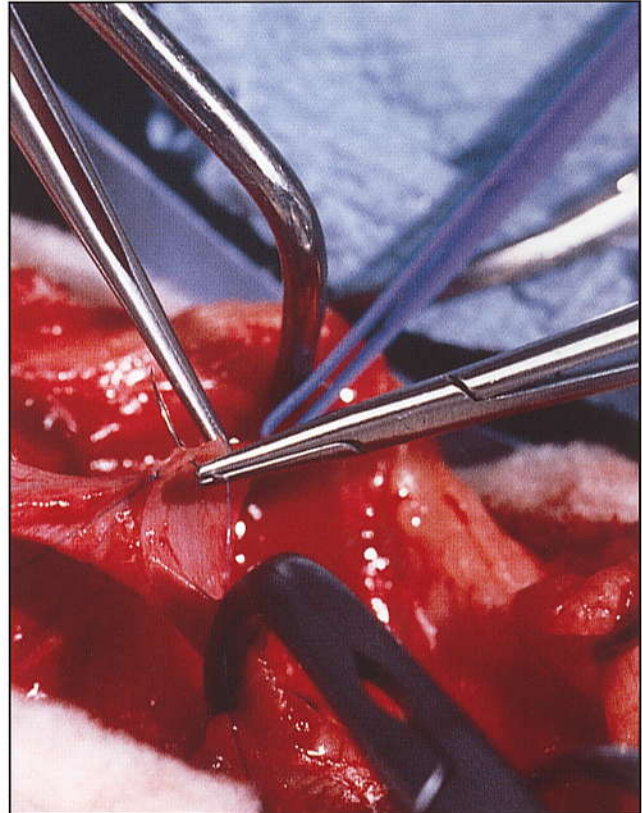

*Custom Surg-I-Loop® silicone loops sizes available upon request.*

### CATALOG NUMBERS:

(2 loops per sterile package, 10 packages per box)

1001-75	red, mini	1001-79	yellow, mini
1001-76	red, maxi	1001-80	yellow, maxi
1001-77	blue, mini	1001-81	white, mini
1001-78	blue, maxi	1001-81M	white, micro
1001-78A	blue, super maxi	1001-82	white, maxi

**Length:** 16", 40 cm

**Width:** (See below)

Super Maxi 5 x 1.27 mm, 3/16" (Blue)

Maxi 2.5 x 1.02 mm, 1/8" (Yellow)

Mini 1.5 x .88 mm, 3/32" (Red)

Micro .5 x .26 mm, 1/16" (White)



For more information or to place an order, please contact your local Scanlan representative, or call...

U.S. & Canada: 800-328-9458 • International: 651-298-0997 • FAX: 651-298-0018 • [www.scanlaninternational.com](http://www.scanlaninternational.com)

**SCANLAN** ★

INTERNATIONAL

*Surgical products designed and manufactured by the Scanlan family since 1921*  
ONE SCANLAN PLAZA • SAINT PAUL, MINNESOTA 55107 • U.S.A.

MK-B-036 Rev. C  
© 2004 Scanlan International, Inc.